

Cite this: DOI: 10.1039/xxxxxxxxxx

# A Method for Studying the Diffusion of Quaternary Ammonium Cations Through Polyelectrolyte Phases<sup>†</sup>

Alexander M. Papiez,<sup>\*a</sup> Justin J. B. Perry,<sup>b‡</sup> and Les R. Dix<sup>a</sup>

Received Date

Accepted Date

DOI: 10.1039/xxxxxxxxxx

www.rsc.org/journalname

The mobility of organic cations in polymeric phases is an important property to consider when using these materials as active ingredients in coatings. Here we describe a method for extracting such compounds from polymeric samples and how analysis of these extracts can yield insights about the diffusivity of molecules in a polymeric phase.

## 1 Introduction

Diffusion in polymeric phases is an important phenomenon which influences many fields. The ability to control the release of active compounds from a polymeric vehicle may be influenced by the diffusivity of these compounds, especially when strong interactions exist between the active compound and the vehicle. Particularly interesting are those cases in which these interactions can be modified to tune the diffusivity of the mobile compound.

In order to assess how different structural features assess the diffusivity of an analyte, the kinetics of analyte release must be measured. The methods used to effect this measurement are highly dependant on both the nature and quantity of the analyte of interest.


**Table 1** Some typical methods used to detect different types of analyte

Analyte	Detection Method	Sensitivity
Transition metal	Flame Photometry	10-1000 ppm
	Flame AAS	1-100 ppm
	Flame AES	< 1 ppm
Organic Cations	HPLC-MS	10ppb - 10 ppm
	GC-MS	10ppb - 10 ppm
	qNMR	10 - 1000 ppm

While HPLC-MS and GC-MS are by far the most sensitive techniques for the investigation of organic compounds, much time and effort must be spent developing and optimizing analyte extraction, pre-concentration and detection methods. Alternative methods, such as flame-photometry and quantitative NMR spectroscopy require less method development but are concomitantly

less sensitive. The aim of this work was to establish whether the kinetics of the ion exchange of quaternary ammonium cations could be studied using a combination of flame photometry and qNMR.

### 1.1 Ion Exchange

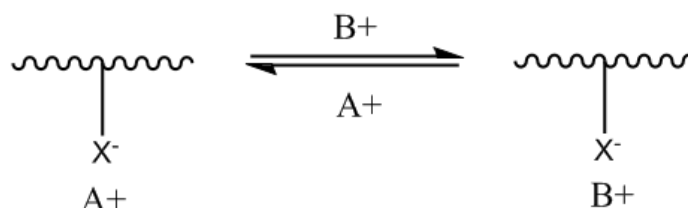

**Fig. 1** Basic outline of the ion exchange equilibrium, where an incoming ion B+ can displace a resin associated ion A+

### 1.2 Basic Principle

Ion exchange describes the phenomenon in which an ion exchange material in some initial ionic form will, in contact with a solution containing some ion of a different type to that already contained within the resin, will take-up that ion while releasing the initial ion.



This phenomenon is of great utility in a large number of industrially important processes; notably the preferential extraction of radioactive isotopes from the waste produced by nuclear reactors.

### 1.3 Ion Exchange materials

A typical ion exchange material consists of a crosslinked polymeric matrix containing acidic or basic side chains (depending upon whether cation or anion exchange is the desired be-


<sup>a</sup> Address, Address, Town, Country. Fax: XX XXXX XXXX; Tel: XX XXXX XXXX; E-mail: xxx@aaa.bbb.ccc

<sup>b</sup> Address, Address, Town, Country.

<sup>†</sup> Electronic Supplementary Information (ESI) available: [details of any supplementary information available should be included here]. See DOI: 10.1039/b000000x/

<sup>‡</sup> Additional footnotes to the title and authors can be included e.g. 'Present address:' or 'These authors contributed equally to this work' as above using the symbols: ‡, §, and ¶. Please place the appropriate symbol next to the author's name and include a \footnotetext entry in the the correct place in the list.

haviour of the resin). One of the most popular co-polymers used as a matrix for ion-exchange material is styrene-divinylbenzene (SDVB)(see figure 2).

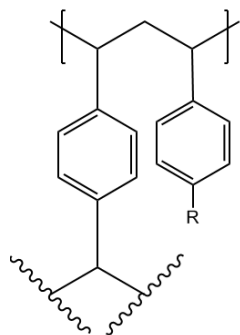


**Fig. 2** The styrene-divinylbenzene co-polymer is one of the most popular skeletons for ion exchange resins. R may be any acidic or basic group, dependant upon whether cation or anion exchange functionality is desired

SDVB based resins are particularly attractive due to the ease with which polymer beads of a well controlled size distribution may be obtained. This is achieved by inverse phase suspension free radical polymerization of the monomers in water. As we shall later discuss, this size control is critical for the production of ion exchange resins which behave in a well defined and predictable manner.

#### 1.4 The Kinetics of Ion Exchange

The kinetics of ion exchange are well understood, with the first pioneering studies undertaken by Boyd at the beginning of the 20th century<sup>1</sup>. The early consensus which was established is that diffusion of ions is the rate-controlling step in ion exchange reactions. There are however two separate diffusive mechanisms which may dominate, depending on reaction conditions, the structure of the ion exchange resin and the nature of the ions undergoing exchange.

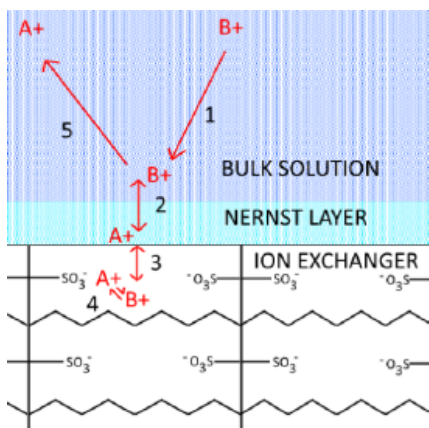

**Fig. 3** Ion exchange consists of a number of steps; 1+5: diffusion of ions through the bulk solution, 2: diffusion of ions through the Nernst Layer, 3: diffusion of ions through the polymeric phase, 4: exchange of ions at the ionic polymer side chain

Figure 3 shows the various sub-processes which make up the overall ion exchange process. It is widely accepted that of all these processes, only diffusion through the Nernst layer and the ion exchanger itself are rate controlling (processes 2 and 3 respectively).<sup>2</sup>

When studying the ion exchange reaction it is important to determine which diffusion type is rate controlling; not only are they subject to significantly different mathematical interpretations, they shed light on diffusion in different regions. In other words, to study diffusivity of a molecule within a polymeric phase, particle diffusion not film diffusion must be rate controlling. Table 2 shows how the reaction conditions can be modified in order to influence which process will be rate controlling.

**Table 2** The effect of different reaction conditions on the two potential rate controlling diffusion types

Condition	Particle Diffusion	Film Diffusion Control
Ion mobility in particle	$\propto$ mobility	No effect
In bulk solution	No effect	$\propto$ mobility
Particle size	$\frac{1}{r^2}$	$\propto \frac{1}{r}$
Capacity of exchanger	no effect	$\propto \frac{1}{X}$
Nature of ionic groups	$\propto$ strength of electrostatic interaction	No effect
Degree of cross-linking	$\propto \frac{1}{\text{crosslinkdegree}}$	No effect

#### 1.5 Measuring Ion Exchange Kinetics

Boyd proposed the following set of equations for the ion exchange reactions (table 3).

**Table 3** Equations describing the ion exchange reaction as stated by Boyd<sup>1</sup>.  $D^i$  = Inter-diffusion coefficient,  $F$  = Fractional approach to equilibrium,  $r_0$  = ion exchange particle radius,  $\Delta r_0$  = thickness of the Nernst layer,  $t$  = time.

Description	Model
Film diffusion in a sphere	$F = 1 - \frac{6}{\pi^2} \sum_{n=1}^{\infty} \frac{1}{n^2} e\left(-\frac{D^i \pi^2 n^2 t}{r_0^2}\right)$
Film diffusion in a polymer slab	$F = 1 - \frac{8}{\pi^2} \sum_{n=1}^{\infty} \frac{1}{(2n-1)^2} e\left[-\frac{(2n-1)^2 \pi^2 D^i t}{4\Delta r_0^2}\right]$
Diffusion through the Nernst Layer	$\log(1-F) = -\left(\frac{[\frac{3D^i}{2.303}]\right)t$

Thus measurement of the extent to which ion exchange reactions have reached equilibrium will permit estimation of the inter-diffusion coefficients for the ion exchange reactions on the involved resin systems. The fractional approach to equilibrium can be monitored by following the ionic composition of the bulk solution or the ion exchange resin; the advantages and disadvantages of each approach are summed up in table 4

**Table 4** Both strategies for studying ion exchange reactions have benefits and drawbacks

Strategy	Positive Aspects	Negative Aspects
Bulk electrolyte monitoring	Compatible with industrial practice Very precise and accurate Easily fits required constraints	Requires complicated experiment Ambiguity about effect of column flow dynamics May require sample preconcentration
Particle content monitoring	Less complicated instrumentation No sample preconcentration required	Requires many distinct observations Incompatible with industrial practice

Where possible, the best approach is of course to combine particle composition analysis with electrolyte composition analysis to come up with a comprehensive trace of the concentrations of the various analytes throughout the reaction. However where this is not possible, one can still deduce the concentration profile of the unmeasured as the velocity of ingress and egress are identical.

To study the ion exchange process a number of steps must be carried out:

- Controlled exposure of a sample of ion exchange particles to an electrolyte solution
- Recovery and washing of the bead sample
- Extraction of analytes from the bead sample
- Quantitation of extracted analytes

We endeavoured to adapt or develop techniques for these processes as follows.

## 2 Experimental

### 2.1 Ion Exchange Reactions

The exchange of sodium for hydrogen on a styrene divinylbenzene commercial resin (AG50WX8) will be used as a representative example. 0.200g of  $H^+$  form AG50WX8 were placed in a nylon mesh basket. This basket was attached to the drive shaft of an overhead stirrer by means of a compression fit over a rubber bung located at the end of the shaft. The basket was set to spinning at  $200 \pm 20$  revolution per minute. Into a jacketed reaction vessel was placed 400ml of 0.500M NaCl solution. If the time for immersion was to be shorter than 2 minutes, then the jacketed reaction vessel was manually elevated such that the flange joint was sealed tight and the nylon basket containing the ion exchange resin was completely submerged. If the planned immersion time was to be 2 minutes or greater then a lab jack was used to accomplish same. After the planned immersion time had been reached, the basket was removed from the electrolyte, allowed to drain, then stirred at 200 rpm for 10 seconds in 600ml distilled water in order to remove any supernatant electrolyte on the bead sample. Finally the nylon basket was separated from the rest of the apparatus to effect removal of the beads which were weighed and placed in a glass sample container. This procedure was repeated for each planned immersion period, normally a minimum of 27 immersion times for a given experiment. After every 6th experiment, the electrolyte solution was disposed of and replenished in order to prevent significant loss of diffusible ion concentration.

#### 2.1.1 Interruption Tests

The procedure for an interruption tests is much the same as the procedure for a standard ion exchange reaction. However for one sample, the basket is reimmersed 30 seconds after washing before a second rinse and sample recovery

### 2.2 Analyte extraction

The recovery of sodium from AG50WX8 shall be used as a representative example. A sample of ion exchange particles, already subjected to the ion exchange procedure as described above is prepared for extraction in the following manner. A 10 ml plastic syringe is plugged with a small quantity of glass wool such that the base of the syringe is covered and no solid matter may pass through. Next, the sample is washed into the vial with distilled water to ensure that the entire contents of the vial are transferred. Finally another glass wool plug is added to the top of the syringe such that the sample of ion exchange particles is encased in glass wool.

Next, 20 ml distilled water is passed through this cartridge to ensure any supernatant liquids which might interfere with the analysis are removed. A collection vial, with a volume of no less than 50 ml is placed underneath the cartridge and 20ml 1M HCl followed by 20 ml distilled water is passed through the resin. Finally a moderate vacuum is applied to encourage any remaining liquids to elute from the cartridge.

The resulting eluent is made up to a volume of 50ml volumetrically, whereupon the concentration of analytes of interest are analyzed.

### 2.3 Assessment of $[Na^+]$ by Flame Photometry

#### 2.3.1 1: Preparation of a calibration curve

Some representative samples are taken from the set to be analyzed: Ideally these should be the samples expected to represent the highest and lowest concentrations of analyte in the data set. The response of the flame photometer to these samples, as compared to a 1000 ppm standard solution is used to estimate the appropriate calibration range for the experiment. Typically for AG50WX8 resins this concentration range will be from 0 - 100 ppm  $Na^+$ .

#### 2.3.2 2: Analysis of samples

Samples are aspirated in the following manner: Each sample is aspirated 3 times, following the third measurement of a sample

the highest calibrant and the blank are reaspirated to ensure stability of the scale. All flame photometry work was conducted on a Jenway 2000 Natural gas flame photometer.

## 2.4 Quantitative NMR spectroscopy

All NMR spectra were recorded on a JEOL 400 ECS NMR spectrometer, operating at a Magnetic field strength of 9 Tesla (400 MHz).

### 2.4.1 Preparation of a calibration curve

In all cases, Maleic acid was used as the quantitative NMR standard. Calibration curves for the analysis of tetramethylammonium chloride will be used as a representative example. A solution of 1000 ppm Maleic acid in D<sub>2</sub>O is prepared by dissolving 1 mg maleic acid per ml of D<sub>2</sub>O. A solution of 1000 ppm tetramethylammonium chloride is made up by dissolving 7.4 mg tetramethylammonium chloride per ml of D<sub>2</sub>O volumetrically.

The total volume of the solution in the NMR sample tube is set at 1ml. At this volume 1μl of stock solution per 1000μl corresponds to 1ppm of sample.

Calibrant solutions are made up such that each solution has a maleic acid concentration of 50 ppm, while the concentration of the tetramethylammonium chloride varies from 10-50ppm.

Collection of NMR spectra is conducted with a view to maximizing the signal to noise ratio for the spectrum. Ideally the signal to noise ratio should be 50:1. Signal to noise ratio is improved by increasing the number of scans obtained for a sample.

## 2.5 Photographic particle size analysis

### 2.5.1 Calibration of the scale

The camera is suspended over the high contrast background and the scale card is placed on the backing. A photograph of this scale is taken at the level of focus which will be used for photographic particle samples. From this image a scale can be derived (pixels/mm) allowing determination of particle sizes

### 2.5.2 Particle size analysis

Two cameras were used in the analysis of particle sizes, one USB microscope and one "smart phone" camera. These two devices provided different levels of zoom permitting particles of different sizes to be analyzed. The method of sample analysis was the same regardless of camera. A sample of ion exchange particle is deposited on the High contrast background. Some sample manipulation takes place in order to separate (as far as is possible) beads from one another to permit clean detection of outlines. Several photographs of each sample are taken.

## 3 Results and Discussion

## 4 Conclusions

## References

- 1 G. Boyd, A. Adamson and M. LS, *J. Am. Chem. Soc.*, 1947, **69**, 2836–2848.
- 2 F. G. Helfferich.